



Thinking ahead. Focused on life.





Veraview X800

The New Frontier for Dental Imaging


The Veraview X800 is an all-in-one dental X-ray unit that produces stunning images for panoramic, cephalometric, and CBCT evaluation.

High resolution, this unit offers a minute voxel size of just 80 µm and features an adjustable horizontal X-ray beam for artifact reduction. Two exposure modes offer control and flexibility with a 360° high definition scan, or a faster 180° rotation with reduced dose. With a range of image sizes and unique features, the Veraview X800 has reached the pinnacle of dental imaging technology.






High Resolution Images

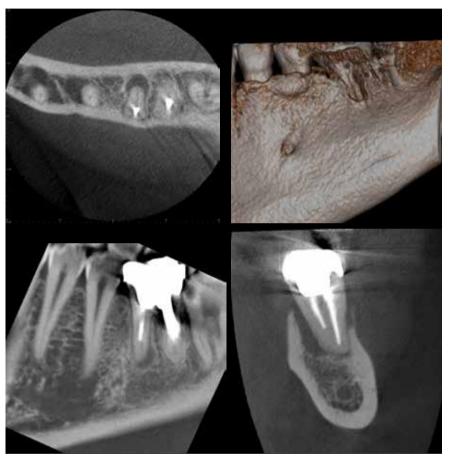



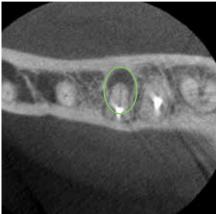

High resolution (80 µm)




Standard resolution (125 µm**)

High-resolution, Limited-field CBCT Imaging

For FOV Ø40 \times H40 exposures, the voxel size is 80 μ m and resolution is 2.5 LP/mm. As shown in the examples above, artifacts are reduced using 80 μ m as compared to 125 μ m voxels.

*Spatial resolution indicates how small objects can be and still be discriminated visually. This is called spatial frequency and is usually expressed as "line pairs per millimeter (LP/mm)". This indicates how many pairs of white and black lines can be discriminated within 1 millimeter; the higher the number, the greater the resolution. MTF (Modulation Transfer Function) is one way to objectively evaluate the line pair resolution and objectively express how many line pairs, and at what level of contrast, can be discriminated. Generally, if MTF is 10%, naked eye discrimination is possible. Spatial resolution does not depend solely on voxel size.


Spatial resolution*

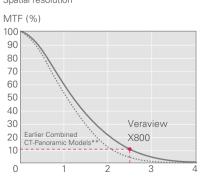

Spatial Frequency (LP/mm)


^{**}Veraviewepocs 3D Series

Horizontal X-ray Beam

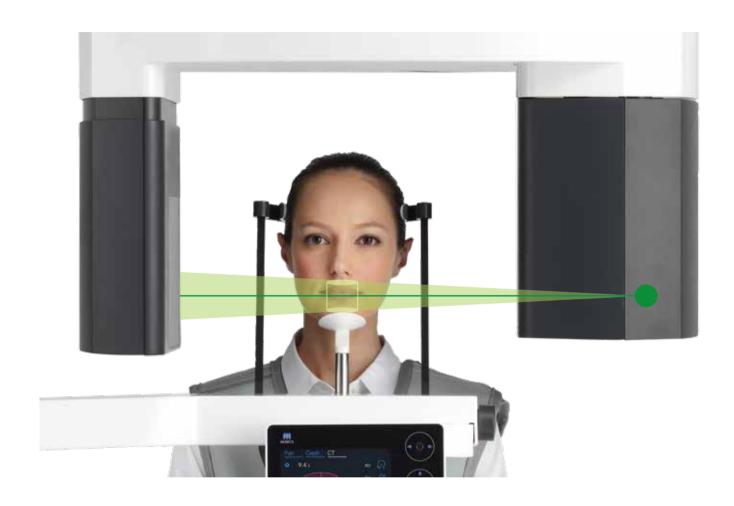


Horizontal X-ray Beam for Minimal Artifacts

For CBCT exposures, the X-ray beam is horizontal during emission which minimizes artifacts and reduces distortion.

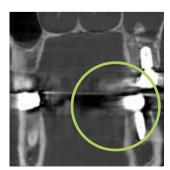

CBCT Image (horizontal beam)

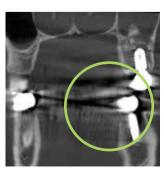


Veraviewepocs 3D series Image (raised beam)


Adjustable for Panoramic Exposures

By shifting the Flat Panel Detector (FPD), the angle of the X-ray beam can be adjusted from horizontal (for CBCT exposures) or raised 5° for panoramic exposures. This slight adjustment suppresses the hard palate during a panoramic and ensures high quality for both image types.



360° Scan

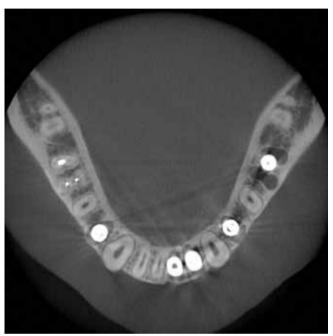

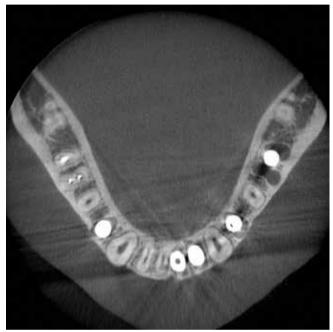


360° Mode Image


180° Mode Image

360° & 180° Exposure Modes

Depending on the diagnostic purpose, the 360° mode can be used for greater detail or the 180° mode can be used for lower X-ray dose and a quicker exposure time.

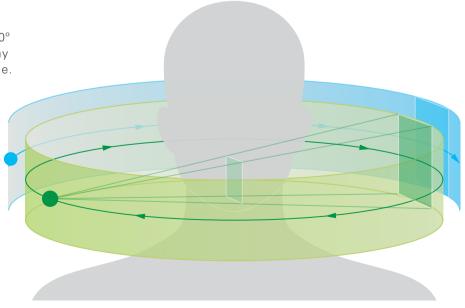

CBCT Positioning with Scout

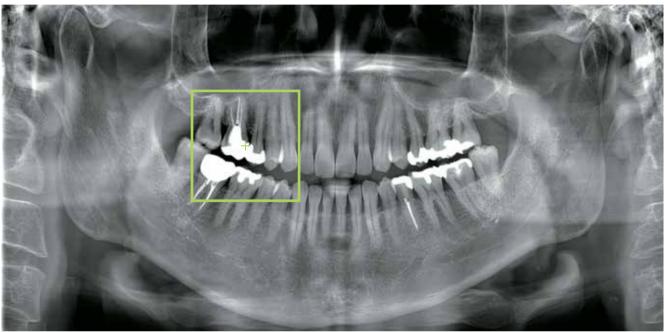


Panoramic Image



Axial Section Image



Sagittal Section Image

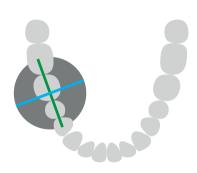




Frontal Section Image

Panoramic Scout

By specifying the region of interest in a panoramic image, positioning and exposure for a limited field CBCT is very simple. This reduces stress for the patient.


After taking a CBCT exposure, double click a cross mark on the panoramic image to display the CBCT data for that region.

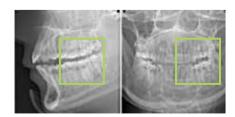


Two-direction Scout

The region of interest is specified by taking lateral and frontal scout images. These images are used to execute accurate positioning for a limited field CBCT exposure.

Various Fields of View

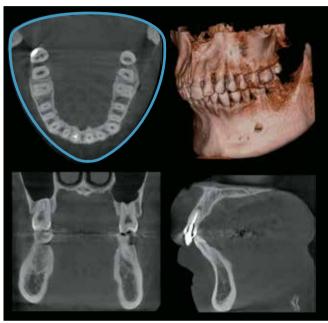


Dental Arch FOV (R 100 x H 80 mm)

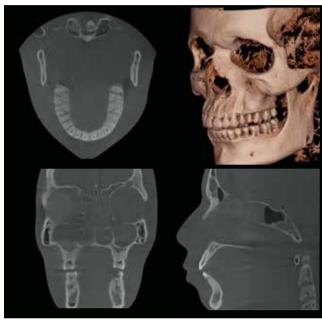

Ø 150 x H 140 mm

Dental Arch FOV Function

A uniquely shaped field of view (FOV) of Ø100 mm encompasses the entire dental arch. Imaging of the entire arch can be executed with less X-ray dose by excluding areas outside the region of interest.

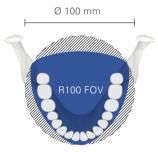

Maximum FOV Ø150

The X800's largest field of view, Ø150, allows a scan of the entire jaw region which is useful for orthodontic, TMJ and occlusal observation and treatment.

Dose Reduction Function

Patient X-ray dose can be reduced by as much as 40%* by lowering the amount of radiation used for areas with greater transparency.

* Compared to when the Dose Reduction function is turned OFF.

Zoom Reconstruction Function**

For the first time, Morita's zoom reconstruction feature is available on a multifunctional unit. After taking an image with a voxel size of 125 µm, reconstruction can be repeated for a higher resolution of 80 µm voxel size without retaking the exposure.

** This function cannot be used for Ø150 exposure.


FOV	Voxel Size	180° Mode	360° Mode	F40	F80*1	R100	F150
\emptyset 40 \times H 40 High Res	0.080 mm	0	0	0	0	0	0
Ø 40 × H 40	0.125 mm						
Ø 40 × H 80							
Ø 80 × H 40	0.125 mm	0	0	_	0	0	0
Ø 80 × H 50							
Ø 80 × H 80							
R 100 × H 40*2	0.125 mm	0	_	_	_	0	0
R 100 × H 50*2							
R 100 × H 80*2							
Ø 150 × H 50*3	0.320 mm	-	0	_		_	0
Ø 150 × H 75*3							
Ø 150 × H 140*4							

^{*1} F80 model is only available in Canada.

^{*2} R100: Dental Arch FOV (Ø100 equivalent)

^{*3} Data used equivalent to 180 degree exposure.

^{*4} Two 360 degree exposures, top and bottom. Data used equivalent to 180 degree exposure.

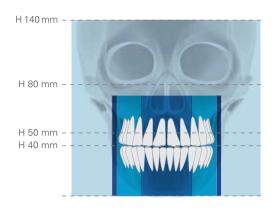



Multiple Fields of View

Veraview X800 offers a total of eleven different FOVs from \emptyset 40 up to \emptyset 150. This unit is appropriate for a variety of specialties and applications such as orthodontics, implantology, periodontics, and endodontics.



Ø 40 x H 80 mm



Ø 80 x H 80 mm



Ø 150 x H 140 mm



Ø 40 x H 40 mm




Ø 80 x H 50 mm



Dental Arch FOV (R 100 x H 80 mm)

Panoramic

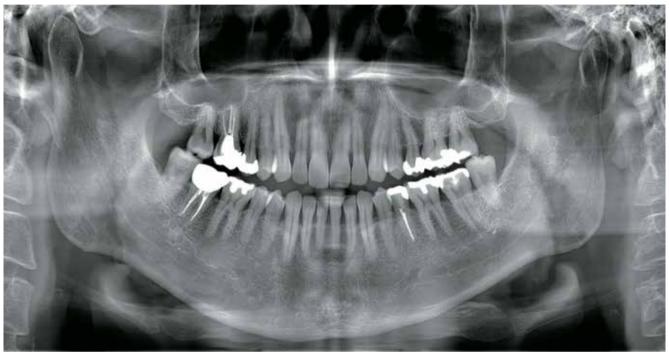


Panoramic Image - Consistent Image Quality and Clarity


AFP + AGS + AIE-HD

Various functions such as AFP (Adaptive Focal Point), AGS (Adaptive Gray Scale), and AIE-HD (Auto Image Enhancement – High Definition) can be combined to obtain images uniformly in focus with the areas of interest clearly observable. Additionally, the clinician may choose between three types of image-layer orbits to suit the individual patient's dentition.

AFP (Adaptive Focal Point)

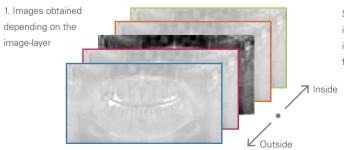


Select the area that is in focus and the entire image will come into focus.

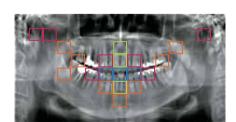

AFP (Adaptive Focal Point)

This function enhances the focus and clarity for the entire scan obtained by the image-layer exposure. Everything in the image from the root apex to the incisor region is in focus.

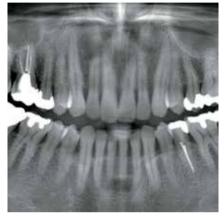

AIE-HD (AFP OFF)



AIE-HD (AFP ON)

AGS (Adaptive Gray Scale)

AGS automatically adjusts density to make the whole panoramic image clearly observable including the dental arch, jaw bone, TMJ, etc.



AFP + AIE-HD (AGS OFF)



AFP + AIE-HD (AGS ON)

AIE-HD (Auto Image Enhancement)

This function optimizes panoramic image processing and makes every detail sharp and clear.




AFP + AGS (AIE-HD OFF)



AFP + AGS (AIE-HD ON)

Panoramic

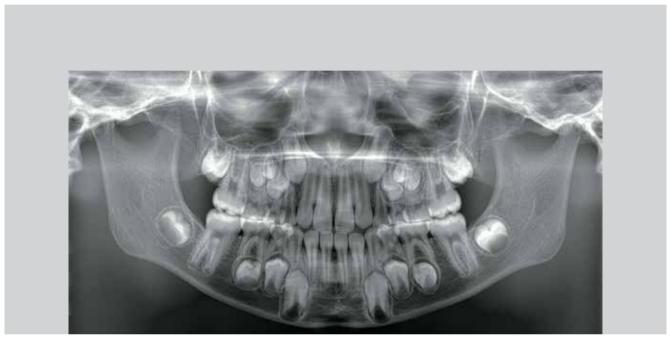


Pediatric Panoramic

For children's smaller jaws, the range of exposure is more narrow to reduce X-ray dose.



Orthographic Panoramic


The X-ray beam intersects the dental arch perpendicularly to reduce overlapping of neighboring teeth.

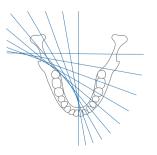


Shadowless Panoramic

This function reduces the shadow caused by the mandibular ramus on the opposite side.



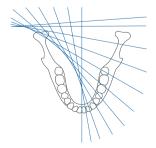

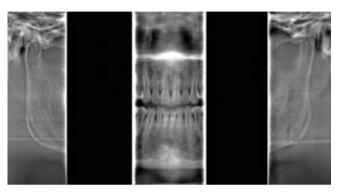

Partial Panoramic


If a complete panoramic image is not needed, a panoramic of only a specified region can be taken. By excluding parts of the dental arch, X-ray dose is reduced.





Bite-wing Exposure

This exposure is useful for prosthetics and diagnosis of mild periodontitis or caries in the proximal spaces of premolars and molars.

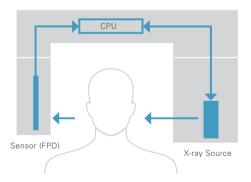

DDAE (Digital Direct Auto Exposure)

During the exposure the flat panel detector detects X-ray transparency in real time and then controls the amount of X-rays emitted to create images with a much better dynamic range.



Image-Layer Form Matches Dental Arch

Three types of image-layer forms are available to suit the individual patient's arch shape.

Cephalometric



High Quality Cephalometric 100 KV tube voltage produces high quality cephalometric images, while high resolution of 96 µm makes it easy to find trace points. Additionally, the soft tissue filter can be adjusted in increments of 5 mm to match the size of the patient.







Partial Cephalometric
Three regions can be
excluded if not needed for
evaluation. This reduces the
patient's X-ray dose.

Face-to-Face Design



Face-to-Face Positioning

Laser beam positioning is more accurate if you have good communication with the patient.

Wheelchair Compatible

The chinrest can be lowered to approximately 34" (865 mm) (short column) to accommodate patients in wheelchairs.

Control Panel

The control panel moves freely so that it can be used from the front, or side, for improved access during patient positioning.



User Interface

The intuitive touch panel is designed for easy operation.










Network System in Clinic

i-Dixel WEB

CBCT and 2D images can be displayed on any PC or tablet computer using a conventional web browser without installing any special software, which is convenient and helpful for patient consultation.

i-Dixel conforms to the following DICOM 3.0

- 1. Modality worklist management service class
- 2. Storage service class
- 3. Modality performed procedure step service class
- 4. Print management service class

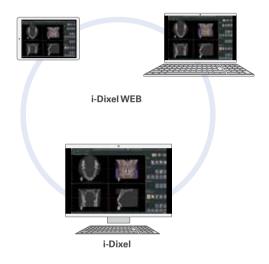







Panoramic / Cephalometric Small Base

Panoramic Large Base

Cephalometric Large Base

Industrial Design: f/p design gmbh

Specifications

Veraview X800 Name:

Model: X800

Order Selection: F40 / F80* / R100 / F150

AC 120V 60Hz Rating:

AC 220 / 230 / 240V 50 / 60Hz

Power Consumption: 2.0 kVA

Weight: Approx. 408 lbs. (185 kg), Approx. 485 lbs. (220 kg)

with cephalometric

Manufacturer: J. MORITA MFG. CORP.

X-ray Tube Voltage: 60 - 100 kV (depending on exposure mode) X-ray Tube Current: 2 – 10 mA (depending on exposure mode)

Nominal Focal Spot: 0.5

Panoramic Exposures: High Speed Mode (standard panoramic) approx. 7.4 sec

High Definition (standard panoramic) approx. 14.8 sec

Panoramic Regions:

Standard Panoramic (Standard, Orthographic, Shadowless), (consistent magnification) Pedodontic Panoramic (Standard, Orthographic, Shadowless),

Maxillary Sinus Panoramic (Anterior, Posterior), Quadruple TMJ, Partial Panoramic, Bite-wing Exposure

Distances measured on a panoramic image are not equal to the

actual distances.

CBCT Exposure Time: Approx. 9.4 sec. (180°) / approx. 17.9 sec. (360°)

CBCT Exposure Regions: F40P / F40CP

- Ø40 x H40, Ø40 x H80

F80P* / F80CP*

- Ø40 x H40, Ø40 x H80

- Ø80 x H40, Ø80 x H50, Ø80 x H80

R100P / R100CP

- Ø40 x H40, Ø40 x H80

- Ø80 x H40, Ø80 x H50, Ø80 x H80

- R100 x H40, R100 x H50, R100 x H80

F150P / F150CP

- Ø40 x H40, Ø40 x H80

- Ø80 x H40, Ø80 x H50, Ø80 x H80

- R100 x H40, R100 x H50, R100 x H80


- Ø150 x H50, Ø150 x H75, Ø150 x H140


Cephalometric

F40CP / F80CP* / F100CP / F150CP Exposures: LA 220 x 250, PA 220 x 200 mm Direction and Size:


> Wear protective aprons and coverings as necessary during X-ray exposure. The unit must be fixed to the floor and wall when installed. If minimal layout dimensions are used, there may be very little space to move around inside


the booth

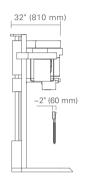

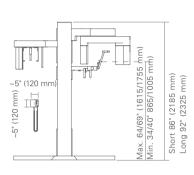



X800 Pan / Ceph / CBCT Units

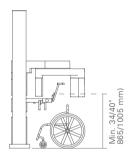

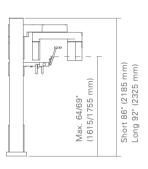

X800 Pan / CBCT Units

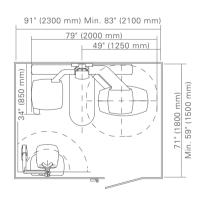

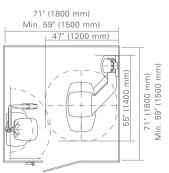

Diagnostic and Imaging Equipment

Treatment Units

Handpieces and Instruments

Endodontic System

Laser Equipment

Laboratory Devices

_

Educational and Training Systems

Auxiliaries



Development and Manufacturing

J. MORITA MFG. CORP.

680 Higashihama Minami-cho, Fushimi-ku, Kyoto 612-8533, Japan T +81. (0)75. 611 2141, F +81. (0)75. 622 4595

Morita Global Website www.morita.com

Distribution

J. MORITA CORP.

3-33-18 Tarumi-cho, Suita-shi, Osaka 564-8650, Japan T +81. (0)6. 6380 1521, F +81. (0)6. 6380 0585

J. MORITA USA, INC.

9 Mason, Irvine CA 92618, USA T +1. 949. 581 9600, F +1. 949. 581 8811

J. MORITA EUROPE GMBH

Justus-von-Liebig-Strasse 27a, 63128 Dietzenbach, Germany T +49. (0)6074. 836 0, F +49. (0)6074. 836 299

MORITA DENTAL ASIA PTE. LTD.

3 Science Park Drive, #01-05 The Franklin Singapore Science Park1, Singapore 118223 T +65. 6779. 4795, F +65. 6777. 2279

J. MORITA CORP. AUSTRALIA & NEW ZEALAND

Suite 2.05, 247 Coward Street, Mascot NSW 2020, Australia T +61. (0)2. 9667 3555, F +61. (0)2. 9667 3577

J. MORITA CORP. MIDDLE EAST

4 Tag Al Roasaa, Apartment 902, Saba Pacha 21311 Alexandria, Egypt T +20. (0)3. 58 222 94, F +20. (0)3. 58 222 96

J. MORITA CORP. INDIA

Filix Office No.908, L.B.S. Marg, Opp. Asian Paints, Bhandup (West), Mumbai 400078, India T \pm 91-22-2595-3482

J. MORITA MFG. CORP. INDONESIA

28F, DBS Bank Tower, Jl. Prof. Dr. Satrio Kav. 3-5, Jakarta 12940, Indonesia T +62-21-2988-8332, F +62-21-2988-8201

SIAMDENT CO., LTD.

71/10 Mu 5, Thakham, Bangpakong, Chachuengsao 24130, Thailand T +66. 38. 573042, F +66. 38. 573043 www.siamdent.com

Catalog Design: f/p design, Germany